


Triathlon®

Trial Cutting Guide (TCG)

Design rationale

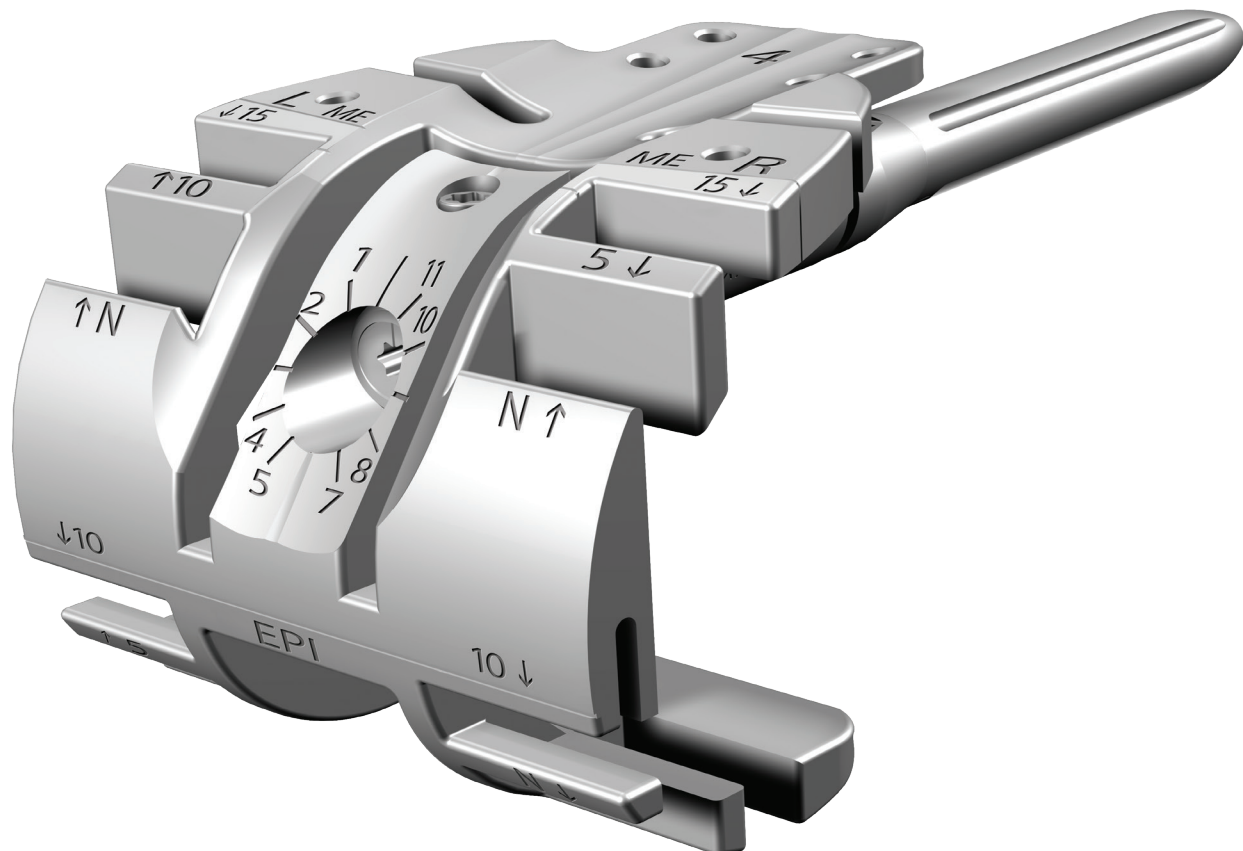

Triathlon Trial Cutting Guide

Component positioning. Simplified.

The Triathlon Trial Cutting Guide (TCG) is designed to simplify femoral bone preparation for a Triathlon Revision knee.¹


The positioning of the femoral component may impact the functional outcomes in primary and revision knee arthroplasty.²⁻¹¹

With the TCG, a surgeon may fine-tune the joint line, gap balancing and femoral rotation—all before committing to any bone cuts.¹

After trial reduction and evaluation through the range of motion, all cuts may be performed and stability can be evaluated with Triathlon PS or TS trials.

Joint line. Where you want it.

Accurate location of the joint line is critical in TKA.^{6,7} Simply align the **medial epicondyle** with the ME scribe mark on the medial side of the TCG and pin through the open slot. This allows all bone cuts to be referenced from the ME, which may assist in accurately and reproducibly relocating the joint line.⁷ Use the cutting slots as a reference to prepare for distal augments, which have been shown to enhance the ability to restore the joint line.⁸

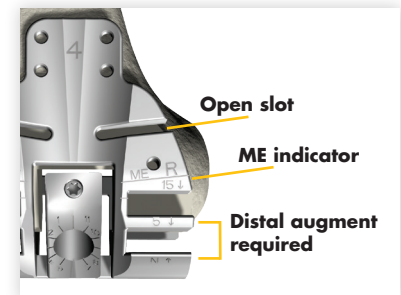


Femoral rotation. Dialed in.

Femoral component rotation impacts patella tracking, and malrotation may result in pain⁹ and complications.¹⁰ Evaluate rotation by performing a trial reduction with the TCG **before committing** to any cuts. Use EPI indicators in extension and flexion to confirm alignment. **Visualize** and perform augment cuts through the instrument; the use of posterior augments may improve accuracy of femoral rotation.¹¹



Gap balancing. Dynamic.

Instability is a leading cause of revisions.¹² Restore flexion and extension gaps with **dynamic offsetting**, which may enhance stability and range of motion in revision TKA.¹³ Offset the component before making augment cuts to posteriorize the component as necessary. **Restoring posterior condylar offset** has been shown to improve flexion stability and be an independent predictor of positive outcomes.¹⁴

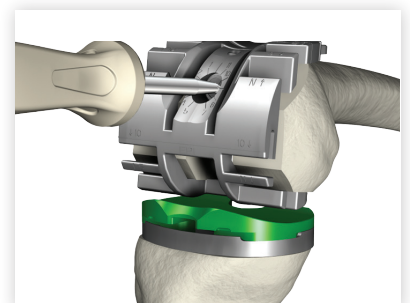


Anatomic component placement

With the Triathlon TCG:

Component positioning

Anterior runout allows you to visualize the final implant positioning

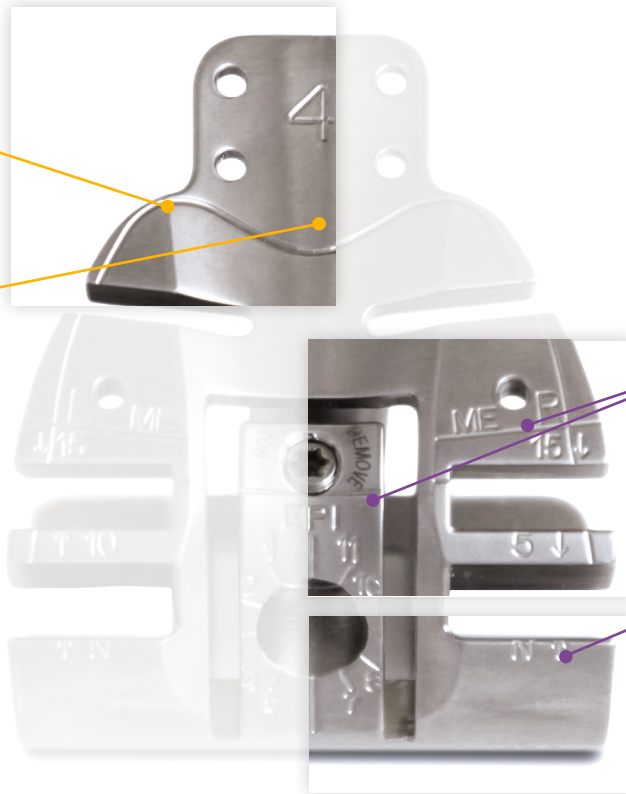

Component positioning

Evaluate patella tracking by referencing the runout of the trochlear groove

Joint line

Locate the joint line 28 mm from the ME by referencing the scribe mark

Joint line

Identify augments needed to restore the joint line, posterior condylar offset and femoral rotation

Gap balancing


Dynamically offset the component to establish flexion gap

Femoral rotation

Locate the ME and transepicondylar axis

Femoral rotation

Identify augments needed to restore rotation/PCO

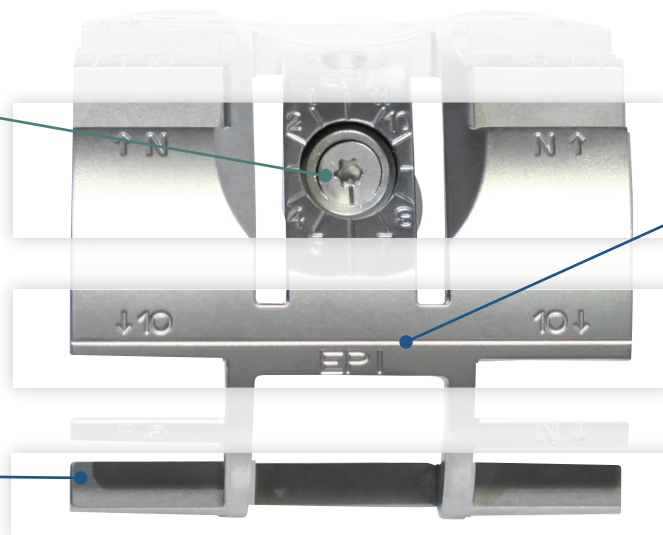

Revision. Simplified.

With the TCG, all femoral bone preparation is performed with one instrument. Specially designed TCG trial inserts provide the ability to assess positioning and stability without insert constraint. Trial femoral components and trial PS or TS inserts can then be used to assess final placement and constraint.¹

In the presence of metaphyseal bone loss, Triathlon Cones are designed to provide a platform for metaphyseal fixation while maintaining desired alignment of the revision components.¹⁵

The Triathlon Cones reaming prep is designed to be simple and quick and allows for an intimate, line-to-line fit and accurate cone positioning.¹⁵⁻¹⁹ Because bone is machined to match the cone, the Triathlon Cones have demonstrated minimized micromotion compared to traditional cones.²⁰

Reduced cone cross sections permit use of larger diameter stems and freedom in regards to implant placement.¹⁵



Revision knees are challenging. From femur prep and bone loss management to PS and TS constraint options, we help simplify your revisions.

References:

1. Triathlon Revision Knee protocol. TRITS-SP-1; 2015.
2. Cummings JF, Carpenter CW, Grood ES, Leach DU, Manley MT. Joint line elevation of a total knee replacement results in reduction of knee flexion. Presented at: Orthopaedic Research Society (ORS) 36th Annual Meeting; February 5-8, 1990; New Orleans, LA.
3. Partington PF, Sawhney J, Rorabeck CH, Barrack RL, Moore J. Joint line restoration after revision total knee arthroplasty. *Clin Orthop Relat Res.* 1999;Oct(367):165-171.
4. Lyons ST, Hofmann AA, Camargo M, Moen C, Feign M. Restoration of the joint line based on the distal femur in revision total knee arthroplasty. Presented at: American Academy of Orthopaedic Surgeons Annual Meeting; March 15-18, 2000; Orlando, FL.
5. Pietsch M, Hofmann S. Early revision for isolated internal malrotation of the femoral component in total knee arthroplasty. *Knee Surg Sports Traumatol Arthrosc.* 2012;20(6):1057-1063. doi:10.1007/s00167-011-1637-3
6. Masini MA, Kester MA. The joint reduction method of revision total knee arthroplasty. *Orthopedics.* 2004;27(8):813-816.
7. Mason M, Belisle A, Bonutti P, Kolisek FR, Malkani A, Masini M. An accurate and reproducible method for locating the joint line during a revision total knee arthroplasty. *J Arthroplasty.* 2006;21(8):1147-1153. doi:10.1016/j.arth.2005.08.028
8. Porteous AJ, Hassaballa MA, Newman JH. Does the joint line matter in revision total knee replacement? *J Bone Joint Surg Br.* 2008;90(7):879-884. doi:10.1302/0301-620X.90B7.20566
9. Barrack RL, Schrader T, Bertot AJ, Wolfe MW, Myers L. Component rotation and anterior knee pain after total knee arthroplasty. *Clin Orthop Relat Res.* 2001;Nov(392):46-55. doi:10.1097/00003086-200111000-00006
10. Berger RA, Cosssett LS, Jacobs JJ, Rubash HE. Malrotation causing patellofemoral complications after total knee arthroplasty. *Clin Orthop Relat Res.* 1998;Nov(356):144-153. doi:10.1097/00003086-199811000-00021
11. Anderson L, Mason JB. Modular augmentation in revision total knee arthroplasty. In: Bono JV, Scott RD, eds. *Revision Total Knee Arthroplasty.* Springer, Cham; 2018:117-130. doi:10.1007/978-3-319-67344-8_9
12. Namba RS, Cafri G, Khatod M, Inacio MCS, Brox TW, Paxton EW. Risk factors for total knee arthroplasty aseptic revision. *J Arthroplasty.* 2013;28(8 Suppl):122-127. doi:10.1016/j.arth.2013.04.050
13. Mahoney OM, Kinsey TL. Modular femoral offset stems facilitate joint line restoration in revision knee arthroplasty. *Clin Orthop Relat Res.* 2006;May(446):93-98. doi:10.1097/01.blo.0000214425.44582.6b
14. Clement ND, MacDonald DJ, Hamilton DF, Burnett R. Posterior condylar offset is an independent predictor of functional outcome after revision total knee arthroplasty. *Bone Joint Res.* 2017;6(3):172-178. doi:10.1302/2046-3758.63.BJR-2015-0021.R1
15. Triathlon Revision surgical protocol. TRITS-SP-2_Rev-2_10451; 2016.
16. Femoral bone prep tolerance analysis. Project 195725:document A0004384; 2015.
17. Tibial bone prep tolerance analysis. Project 195725:document A0004385; 2015.
18. Stryker technical report: Report number RD-14-050. Micromotion of Triathlon Tibial Cone augment constructs. 2015.
19. Stryker technical report: Report number RD-14-049. Triathlon Tritanium Femoral Cone augment micromotion. 2015.
20. Alipit VL, Kirk AE, Stoker M, Meneghini RM. Quantitative stability assessment of novel porous metal metaphyseal femoral cones for revision knee arthroplasty. Presented at: Orthopaedic Research Society (ORS) Annual Meeting; March 5-8, 2016; Orlando, FL.

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented is intended to demonstrate the breadth of Stryker's product offerings. A surgeon must always refer to the package insert, product label and/or instructions for use before using any of Stryker's products. The products depicted are CE marked according to the Medical Device Regulation 2017/745 or the Medical Device Directive 93/42/EEC. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your sales representative if you have questions about the availability of products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Stryker, Triathlon. All other trademarks are trademarks of their respective owners or holders.