

Prime Connect™



Smart stretcher. Confident care.

Empower your team with Prime Connect—Stryker’s first smart, connected stretcher. With both visual and audible bed exit alarms, Prime Connect can communicate with the nurse call system to help prevent patient injury from falls.




Enhanced local monitoring

Prime Connect makes it easy to acknowledge bed exit alarms in loud, fast-paced environments with the ability to illuminate both the corridor light and footboard display when activated.



Bed Exit

For patients who pose a flight or fall risk, bed exit can be set to alarm the nurse when the patient begins to exit the stretcher.

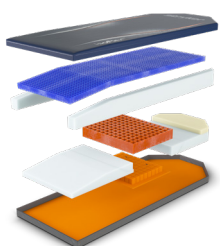

IsoFlex™ SE gel support surface

Offers pressure redistribution and comfort as a component to aid in the treatment and prevention of pressure ulcers.

Advanced mobility options

Big Wheel reduces start-up push force and steering efforts, while Zoom motorized drive nearly eliminates manual pushing.

Recovery chair position

Helps eliminate the need for unnecessary transfers by providing a comfortable seated position.

Lift Assist® backrest

Uses the patient’s weight to raise the head of the stretcher while the drop seat design helps reduce boosting. Also available as a fully electric option.

Patient controls

The patient has control of their own comfort with the staggered control panels on the Glideaway siderails.

Electric-powered lift pedal*

Eliminates manual pumping, allowing for quick patient access.

Power-washable design*

With 318 kg safe working load.

*Models 1105 and 1115 only



Prime Connect

Standard features

- Nurse call integration
- Enhanced footboard
- Bed exit
- Less than 57.7 cm low height (Models 1105 and 1115)
- Litter-mounted oxygen bottle holder
- Center-locking brake system with dual end activators
- Dual, side-mounted foot control hydraulics with uni-lower pedal
- Dual end brake/steer controls
- Glideaway siderails
- Four integrated IV receptacles
- Trendelenburg/reverse Trendelenburg positioning
- 20.3 cm Omni surface Casters with integrated wheel covers
- Pneumatic backrest/stationary foot
- Thermoformed ABS base
- Integrated utility tray
- Pop-up push handles
- Electric powered knee gatch
- Lift Assist backrest design with drop seat
- Three-sided hydraulic controls
- Dual-ended siderail release

Specifications

76.2 cm (30") litter

Model number

5th Wheel Mobility	1105
Big Wheel Advanced Mobility	1115
Zoom Motorized Drive System	1125

Overall length

1105 and 1115	215.9 cm (85")
1125	218.4 cm (86")

Overall width

Siderails up	96.5 cm (38")
Siderails down	77.5 cm (30.5")

Patient surface

76 x 191 cm (30" x 75.25")

Weight capacity

318 kg (700 lb)

Height range (1105, 1115)

High	86.4 cm (34")
Low	57.7 cm (20.75")

Height range (1125)

High	86.4 cm (34")
Low	58.4 cm (23")

Litter positioning

Backrest	0° - 90°
Knee gatch	0 - 13.9 cm (5.5")
Trend./Reverse Trend. (1105, 1115)	+17/-17 (+/- 3)
Trend./Reverse Trend. (1125)	+17/-16 (+/- 3)

Siderail length

36 x 147 cm (14" x 57.75")

Caster diameter

20 cm (8")

EMC

IEC 60601-1-2

Optional features

- IsoFlex SE support surface
- ComfortGel SE support surface
- UltraComfort SE support surface
- Permanent two-stage IV pole
- Serving tray
- Defibrillator tray
- Defibrillator tray/foot extender/charting surface
- Chartholder
- Serving tray holder/footboard
- Siderail pads
- Upright oxygen bottle holder
- Restraint straps

Warranty

2 year parts, labor, travel

Electric functionality

- Electric patient controls
- Foot end operator controls
- Electric lift
- Bed exit
- 120V or 240V available

This document is intended solely for the use of healthcare professionals. A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker or its affiliated entities own, use, or have applied for the following trademarks or service marks: Big Wheel, ComfortGel, Glideaway, IsoFlex, Lift Assist, Prime Connect, Prime Series, Stryker, Vision, Vocera, Vocera Engage, Zoom. All other trademarks are trademarks of their respective owners or holders. The absence of a product, feature, or service name, or logo from this list does not constitute a waiver of Stryker's trademark or other intellectual property rights concerning that name or logo.

