


Sage PrimaFit™ External Urine Management for the Female Anatomy

REF 5400-X

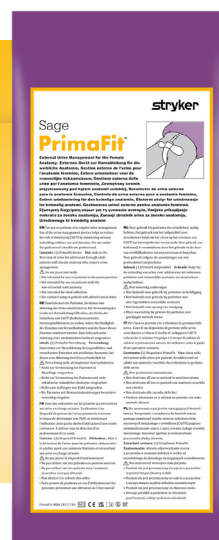

Packaging Information

Contents:

- (1) PrimaFit
- Single use only
- Not made with natural rubber latex
- Non-sterile

Instructions:

Read instructions on package prior to using product

Indications:

For use on patients who require urine management. Use of the urine management devices helps to reduce the risk of developing CAUTI by minimizing urinary indwelling catheter use and duration. For use by healthcare professionals.

Contraindications:

Do not use on patients with urinary retention

Warnings:

Do not place internally

Cautions:

- Not intended for use on patients in the prone position
- Not intended for use on patients with non-retracted male anatomy
- Not intended for use on stool collection
- Use caution if using on patients with altered mental status

Discontinue use:

- When independent or assisted toileting is feasible
- During patient ambulation

Specifications

Product name	Sage PrimaFit External Urine Management for the Female Anatomy
Product code	5400-X
Description	External Urine Management for the Female Anatomy
Manufacturer accreditation	ISO 13485:2016
Expiration date	None
Height of product	Device: 33.0 cm
Length of product	Device: 12.7 cm
Case quantity	40

Sage PrimaFit External Urine Management System for the Female Anatomy

Features and benefits

Sage PrimaFit External Urine Management for the Female Anatomy

Flex-Fit Core:	<ul style="list-style-type: none"> Bends and contours to fit between the labia and against the urethral opening. Maintaining its shape and hold.
Tapered End:	<ul style="list-style-type: none"> Gently fits in the perineum, helping to secure the device in position.
Silicone Adhesive Pad:	<ul style="list-style-type: none"> Silicone based adhesive keeps top of device in place on the patient's supra-pubic area. Release liner tabs allow for quick placement and easy, clean removal from patient's skin.
Ultra-Soft Fabric:	<ul style="list-style-type: none"> Diverts urine away from the skin helping to keep the skin dry.
Latex-free:	<ul style="list-style-type: none"> Reduces the risk of allergic reaction.
Core Structure:	<ul style="list-style-type: none"> Absorbs fluids to be suctioned to collection canister. Diverts fluids away from skin, helping to reduce risk factors associated with IAD.
Tiered Connection Port:	<ul style="list-style-type: none"> Connects device to hospital suction tubing. Allows for easy and secure connection to hospital suction tubing.
Long, flexible tubing:	<ul style="list-style-type: none"> Accommodates a range of body types while allowing multiple suction tubing angles



Sage Products LLC
3909 Three Oaks Road
Cary, IL 60013 USA
800 323 2220 ; +1 815 455 4700

stryker.com

SAGE-UM-OTHD-894194_REV-0_en_us
Copyright © 2024 Stryker
11/2023

This document is intended solely for the use of healthcare professionals.

A healthcare professional must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that healthcare professionals be trained in the use of any particular product before using it.

The information presented is intended to demonstrate the breadth of Stryker product offerings. A healthcare professional must always refer to the package insert, product label and/or instructions for use before using any Stryker product.

Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Stryker, PrimaFit. All other trademarks are trademarks of their respective owners or holders.

The products depicted are CE marked in accordance with applicable EU Regulations and Directives.

PrimaFit is a Class I medical device CE.

This document is not intended for distribution outside EU, UK and EFTA.