


stryker

Pangea™

PeriPRO™ Distal Femur Plate

Design rationale

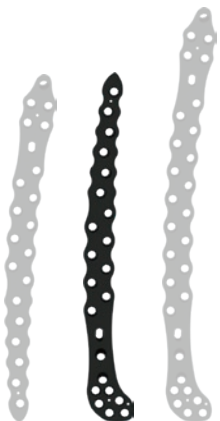

Pangea™ PeriPRO™ Distal Femur Plate

Design rationale

stryker

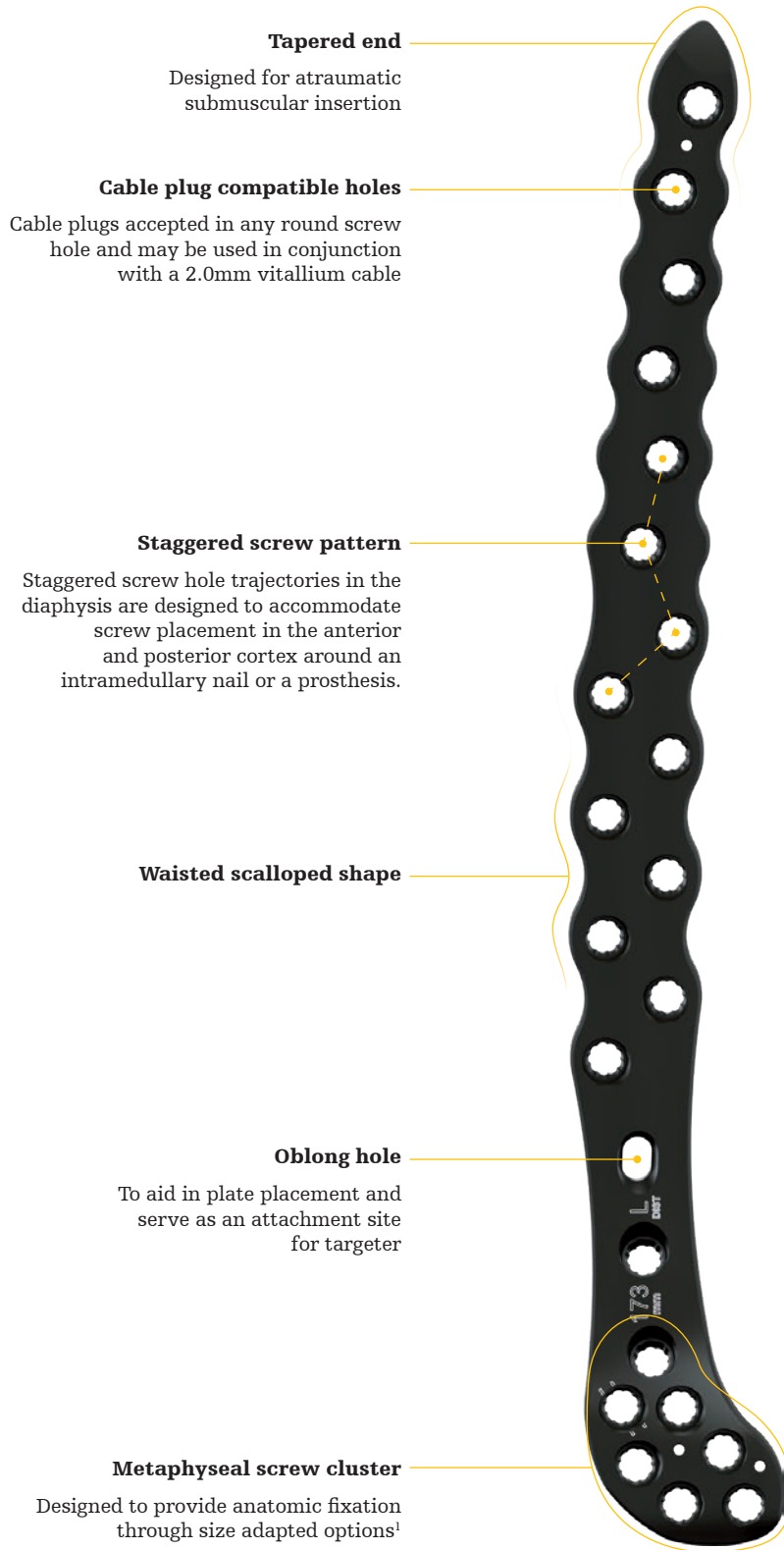


Plate placement

- Radius of Curvature and distal screw cluster plate width changes for longest plate length to facilitate fit¹

Final construct examples



Pangea™ PeriPRO™ Distal Femur Plate

Design rationale

stryker

Plate offering

Length: 173mm

Max. TKA Stem Coverage: N/A



Length: 241mm

Max. TKA Stem Coverage: 100mm



Length: 280mm

Max. TKA Stem Coverage: 120mm



Length: 318mm

Max. TKA Stem Coverage: 160mm



Length: 355mm

Max. TKA Stem Coverage: 200mm



Length: 396mm

Max. TKA Stem Coverage: 240mm




Technical specifications

- Material: Ti6Al4V ELI (plates and non locking screws) and Cobalt chrome (CoCr alloy locking screws)
- Distal femur plate lengths: 173 - 396mm
- Variable angle locking within a 30° cone
- Left and right anatomic plate options

Screw platform

Multiaxial locking	5.0		22-120mm
	4.0		22-95mm
	5.0		10-20mm
Cortex	4.5		22-150mm
Cancellous	6.0		20-150mm
	6.0		30-150mm
	6.0		45-150mm
Cable plug			

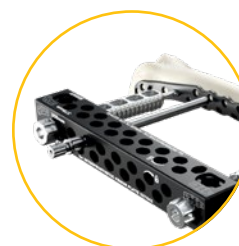

MIS Targeting

Available for all plate types.

Variable angle screw trajectory

The PeriPRO Distal Femur Plates offer an evidence-based design with a variety of plate lengths, screw diameters and SOMA designed screw trajectories, allowing for the ability to achieve a desired construct.¹



References:

1. Internal Report № D0000257417. Rev AA, Selzach, Switzerland

This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery. The information presented is intended to demonstrate the breadth of Stryker's product offerings. A surgeon must always refer to the package insert, product label and/ or instructions for use before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your sales representative if you have questions about the availability of products in your area. Not intended for promotional or marketing use outside of the US and Canada.

Stryker Corporation or its affiliates own, use, or have applied for the following trademarks or service marks: Pangea, PeriPRO, Stryker. All other trademarks are trademarks of their respective owners or holders.

Content ID: PERI-EM-4, 08-2023
Copyright © 2023 Stryker



Manufacturer:

Stryker GmbH
Bohnackerweg 1
2545 Selzach, Switzerland

stryker.com