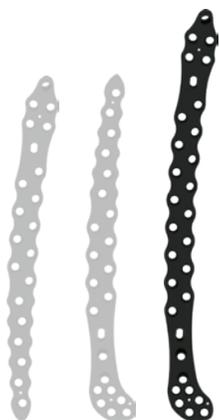
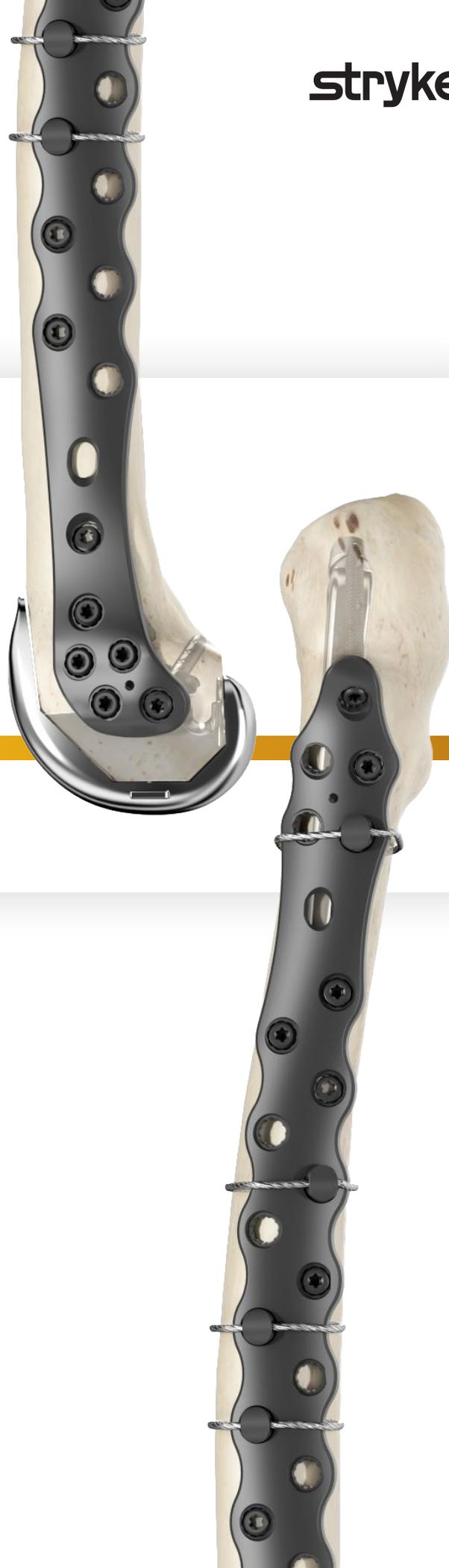


Pangea™

PeriPRO™ Interprosthetic Femur Plate

Design rationale



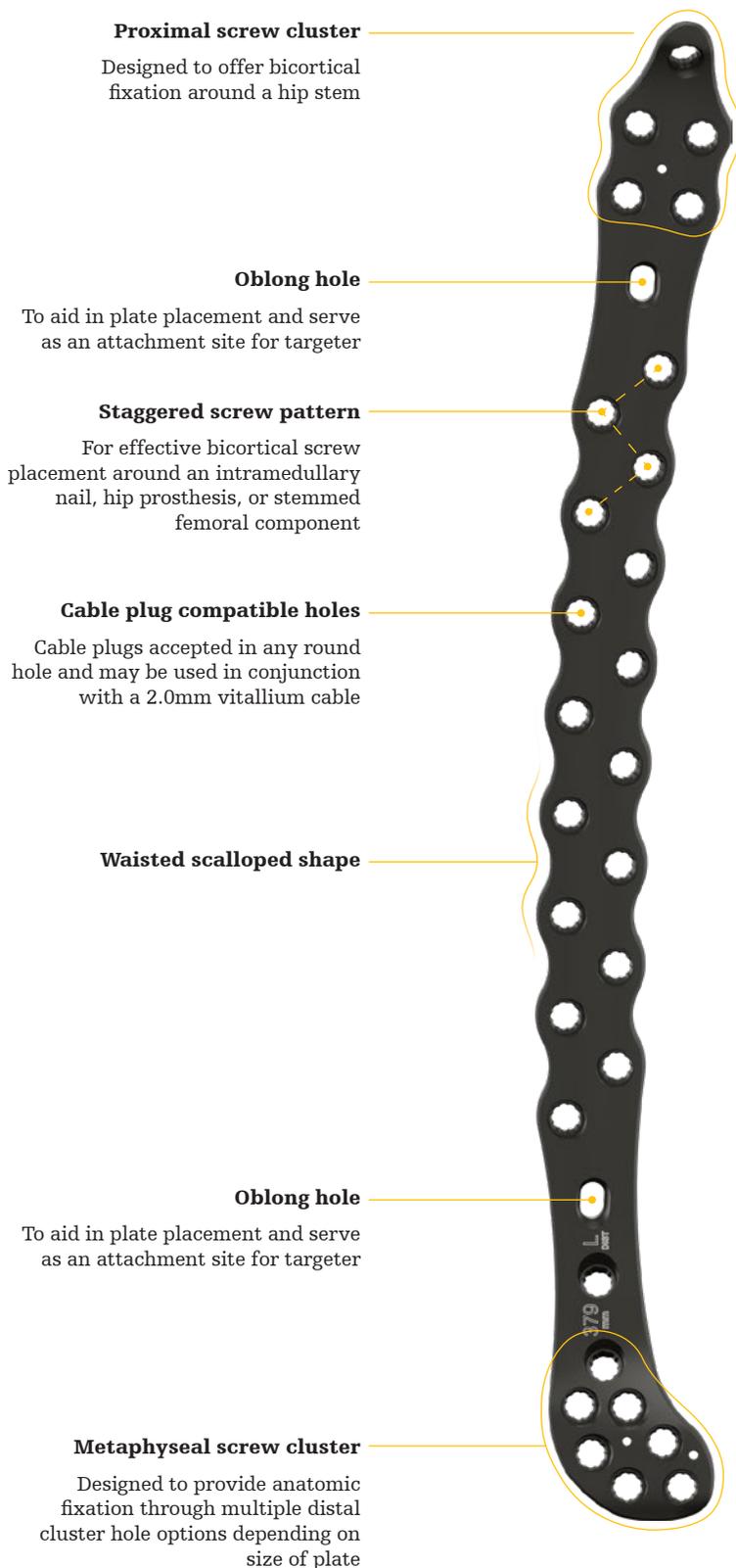


Plate placement

- The first dedicated interprosthetic plate designed to span from the femoral condyle to the greater trochanter to maximize plate working length
- As the length of the plate changes, so does the radius of curvature.

Final construct examples



The first dedicated and patented interprosthetic plate designed to span the entire femur

Pangea™ PeriPRO™ Interprosthetic Femur Plate

Design rationale



Plate offering

Length: 320mm



Length: 334mm



Length: 358mm



Length: 379mm



Length: 402mm



Length: 420mm




Number of holes in the distal cluster varies depending on the size of the plate

Technical specifications

- Material: Ti6Al4V ELI (plates and non locking screws) and Cobalt chrome (CoCr alloy locking screws)
- Interprosthetic femur plate lengths: 320 – 420mm
- Variable angle locking within a 30° cone
- Left and right anatomic plate options

Screw platform

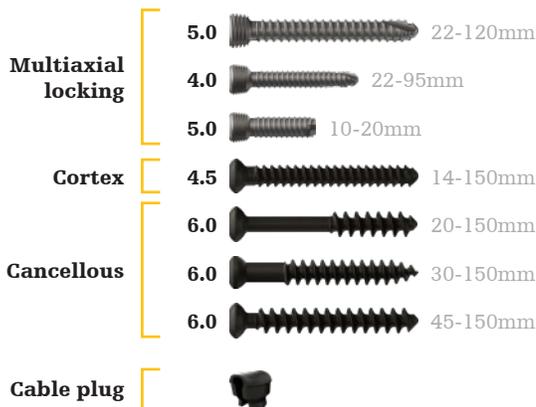

MIS Targeting

Available for all plate types.

Variable angle screw trajectory

Variable angle locking screws are compatible in all circular holes across the entire plate to provide options for bicortical fixation around a prosthesis



This document is intended solely for the use of healthcare professionals. A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery. The information presented is intended to demonstrate the breadth of Stryker's product offerings. A surgeon must always refer to the package insert, product label and/ or instructions for use before using any of Stryker's products. Products may not be available in all markets because product availability is subject to the regulatory and/or medical practices in individual markets. Please contact your sales representative if you have questions about the availability of products in your area. Not intended for promotional or marketing use outside of the United States and Canada.

Stryker Corporation or its affiliates own, use, or have applied for the following trademarks or service marks: Pangea, PeriPro, Stryker. All other trademarks are trademarks of their respective owners or holders.

Content ID: PERI-EM-6, 08-2023

Copyright © 2023 Stryker



Manufacturer:

Stryker GmbH
Bohnackerweg 1
2545 Selzach, Switzerland

stryker.com